


Tap- and Bulk Density Test Instrument - Type PT-TD1

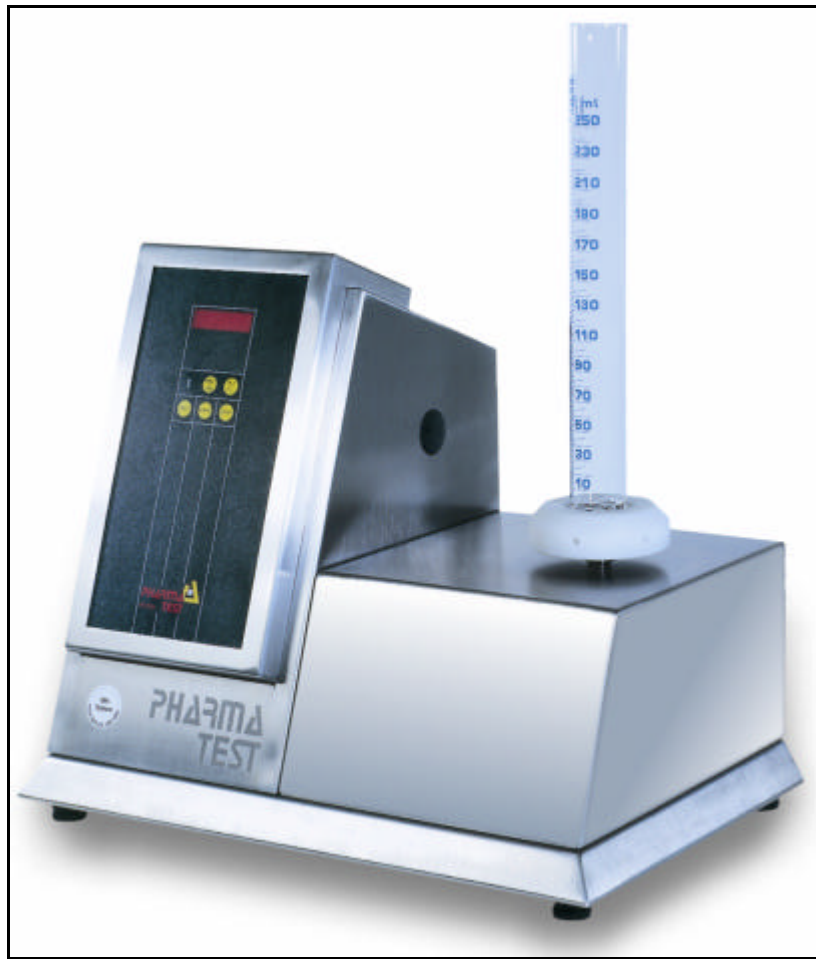

The PT-TD1 testing instrument is designed for tapping or bulk density testing of powders and granules in compliance with USP <616> App. 2, EP <2.9.15> and DIN 53194. The instrument is used to determine additional physical product characteristics with others like particle size, size distribution and shape. Together with information describing the behavioural characteristics such as flow properties, density etc. the results are used to improve the design quality during development as well as the conformity of product in quality assurance and in-process control.

A graduated 250 ml glass cylinder takes the sample to be tested. Driven by a silent DC motor the stroke mechanism performs 250 strokes per minute of a fixed 3 mm height. The total number of strokes or operating time is entered using the keyboard of the apparatus. A digital display shows the elapsed time or number of strokes to the operator who can always recall the nominal information entered previous to the test start.

Usually the sampled is 1250 times tapped, thereafter the volume of the tapped sample is controlled and another 1250 tapping follow, this procedure ist repeated until the difference off the consequent test becomes less than 2%. This result is used to calculate tap volume and tap density. Tap volume is expressed in cm^3/g and calculated as volume divided by weight, tap density is expressed in g/cm^3 and calculated as weight divided by volume.

PHARMA TEST GmbH
Siemensstrasse 5
D-63512 Hainburg (GER)

☎ +49 6182 9532-600
☎ +49 6182 9532-650
✉ email@pharma-test.de
💻 www.pharma-test.com



Options

- 100 ml cylinder and support for less product tests, easy exchangeable with the standard 250 ml set-up (order no. 490-0200)
- Noise Protection Hood



To test for Apparent Volume/Density, or to obtain the Tap Density value for powders and granulates, the sample has to be tapped in a measuring cylinder in order to compact it to a “minimum” volume. The device, by the nature of the test, is relatively loud when in operation. In order to reduce the environmental impact of the noise produced we have introduced a special noise reduction enclosure which reduces the overall sound level by approximately 10dBA. The whole PT-TD 1 fits into the enclosure and is accessed for cleaning and/or sample change through a hinged front door (order no. 49-02100)

Technical Data

Display:	LED Display
Keyboard:	numerical and function keys
Tapping frequency:	250 ± 10 strokes/minute
Tapping height:	3.0 ± 0.1 mm
Number of strokes:	selectable from 1 - 9999
Number of test stations:	1 (to avoid any interference)
Noise Level in operation:	approx. 74 dB (without hood)

Weights and Dimensions

Net weight:	approx. 11 kg
Gross weight:	approx. 20 kg
Packaging:	approx. 60 x 60 x 40 cm

We reserve the right to make technical changes without any prior notice

PHARMA TEST GmbH
Siemensstrasse 5
D-63512 Hainburg (GER)



+49 6182 9532-600
+49 6182 9532-650
email@pharma-test.de
www.pharma-test.com

